

SCIENTEX

Small darkroom

Small type
For microscope

**Make possible to stabilize
of optical measurement**

- For shielding of light on the optical bench
- For the fluorescent microscope
- For safety measure against leak of laser



**For optical measuring
on the desk and table**

Small type B-S-XX

For device that complex form


For microscope B-M-XX

The treatment of preventing dust and static electricity are given to our darkroom.

You can darken your work area locally.

This darkroom can be used in the wide field.

Working with microscope, Chemistry and the field of the material research...etc

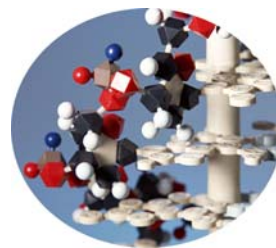

Custom made service

(B-S-XX ~1.5meter)

(B-M-XX ~1meter)

Make up easily by frame pipes!!

Time:10min~30min

(Depends on the size)

Instruction manual

with pictures

Available at clean room

(Only Strong magnet type)

Small type (B-S-XX series)

Available at clean room (B-S-XX-CN)



Hook-and-loop
fastener type

Set up Strong magnets Instead of
Hook-and-loop fastener,



Strong magnet type
(Available at clean room)

For optical measuring on the desk and table

- The edges of curtain are stopped by Hook-and-loop fastener
- For clean room is used strong magnet type

List (Standard)	Available at clean room	Size(mm) *(W)=Entrance
B-S1	B-S1-CN	(W)500 × (D)500 × (H)400
B-S2	B-S2-CN	(W)500 × (D)500 × (H)600
B-S3	B-S3-CN	(W)500 × (D)700 × (H)400
B-S4	B-S4-CN	(W)500 × (D)700 × (H)600
B-S5	B-S5-CN	(W)700 × (D)500 × (H)400
B-S6	B-S6-CN	(W)700 × (D)500 × (H)600
B-S7	B-S7-CN	(W)700 × (D)700 × (H)400
B-S8	B-S8-CN	(W)700 × (D)700 × (H)600
B-S9	B-S9-CN	(W)1000 × (D)1000 × (H)900

We can provide special size on request.
Width: ~1.5meter / Depth & High: ~1meter

Specification

<Frame>diameter 30mm,plastic coating iron pipe

<Weight>6-10kg(the total of the frame and the dark curtain)

<The quality of dark curtain> The shield light ability is the first rank (100-99.99% shielding . Refer to a note 1 for the details.)

<Front side> Polyester 100% Non-flammable material.

<Rear side>Polyurethane black. Laminating processing. The treatment of preventing dust and static electricity.

<Shielding light ability>At Visible ray, Ultraviolet ray and Infrared ray: more than 99.99%.(Refer to a note 1)

*The skirt of dark curtain is made a little long to prevent optical omission.

*Frame pipe and dark curtain are taken apart to ship them. The tools and the manual are packed with them.


For microscope (B-M-XX series)

Available at clean room (B-M-XX-CN)



Edge of dark curtain

Support

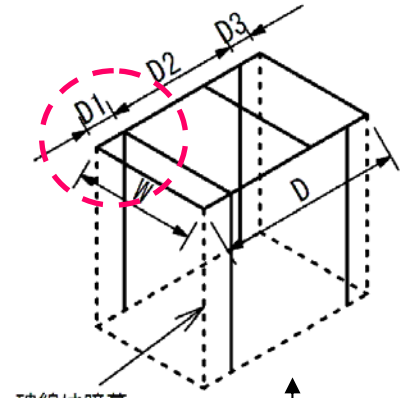

Framework drawing

The space is narrow or For device that complex form

- The edges of curtain are stopped by Hook-and-loop fastener
- For clean room is used strong magnet type
- *The front of dark curtain can be fastened by the Hook-and-loop fastener

List (Standard)	Available at clean room	Size(mm) *(W)=Entrance
B-M1	B-M1-CN	(W)580 X (D)790 X (H)800 (D1) 200, (D2)540, (D3)50 Space to set a dark room (W)580 X (D2)540
B-M2	B-M2-CN	(W)580 X (D)850 X (H)800 (D1) 200, (D2)600, (D3)50 Space to set a dark room (W)580 X (D2)600
<p>We can provide special size on request. Width: ~1.5meter / Depth & High: ~1meter</p>		

*The specification of our product may be changed for the reason such as improvement and the progress of technology without notice.

2015/6/3

SCIENTEX Inc.

Sales department (Business hour 9:00~18:00)
1-4-10-4 Shinmiyakoda Kita-ku Hamamatsu-city
Shizuoka-pref 431-2103 Japan
URL <http://www.scientex.co.jp> E-mail info@scientex.co.jp

Overseas agency